


DuraFlex™ and DuraFlex™ Plus Epidural Catheters



DuraFlex™
catheter

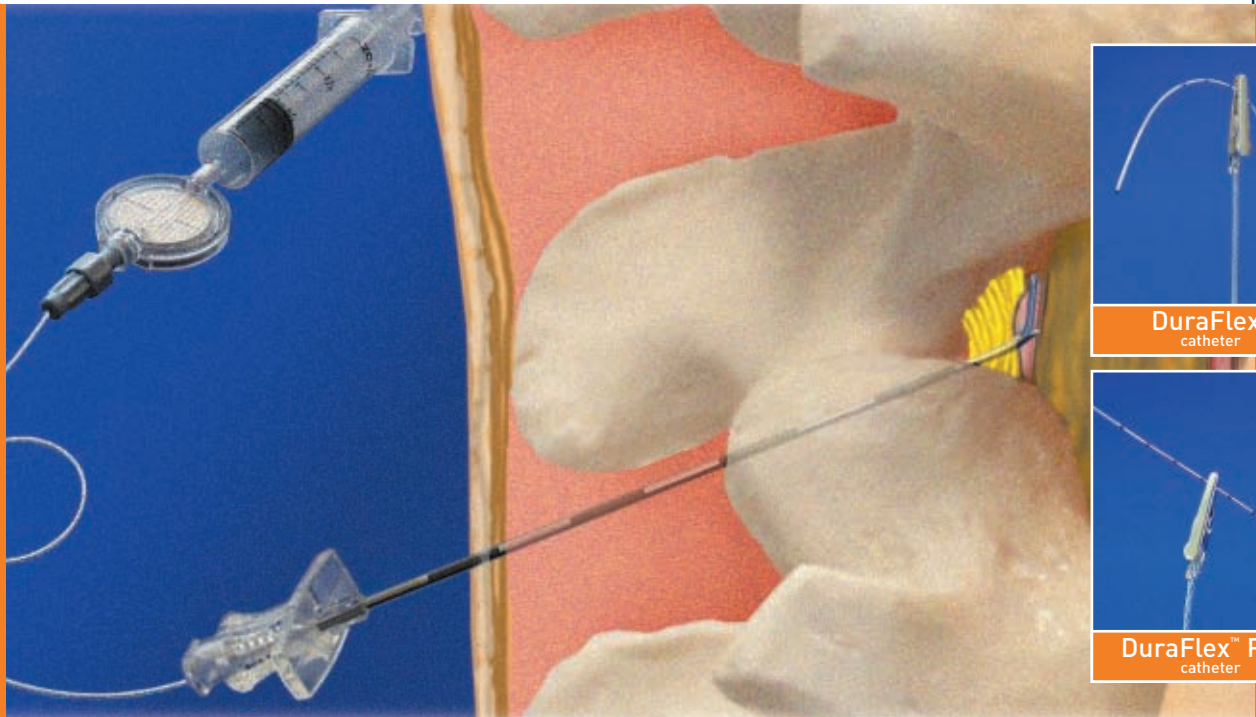

DuraFlex™ Plus
catheter

choices without compromise
in wire reinforced catheters



You now have a choice in epidural wire reinforced catheters from Smiths Medical.

DuraFlex™ Catheters

- > One of the softest catheters on the market

DuraFlex™ Plus Catheters

- > Have the same soft tip as DuraFlex™ catheters, with a stiffer base to prevent coiling and allow easy threading

Both catheter styles combine a flexible wire core enclosed within a soft nylon jacket. The altered spacing of the catheter coils at the distal tip provides increased tip flexibility and enhances the distal tip softness.

Extra flexibility and strength without compromise—where you want it.

DuraFlex™ and DuraFlex™ Plus epidural catheters maintain a balance of strength and flexibility for delivering medication during continuous epidural procedures.

Soft tip.....
easy threading.....
precise stiffness

DuraFlex™ and DuraFlex™ Plus Epidural Catheters


Extra Choice—when you want it.

- > Extra large window opening at two locations on the catheters give better visualization of flash due to inadvertent vascular penetration during catheter placement.
- > A rounded tip instead of a sharp blunt cut softens placement near tender nerves.
- > A secure catheter-to-connector connection prevents leakage. A choice of base stiffness to help ease insertion and placement.
- > Multiport or open-ended tip configuration to ensure medication distribution consistent with therapy requirements.

Extra components—when you want them.

You can choose to incorporate the DuraFlex™ and DuraFlex™ Plus epidural catheters into trays with other components, such as the state-of-the-art Pulsator® L.O.R. syringe.

Your Smiths Medical sales representative will be able to assist you with an optional combination of complementary components.

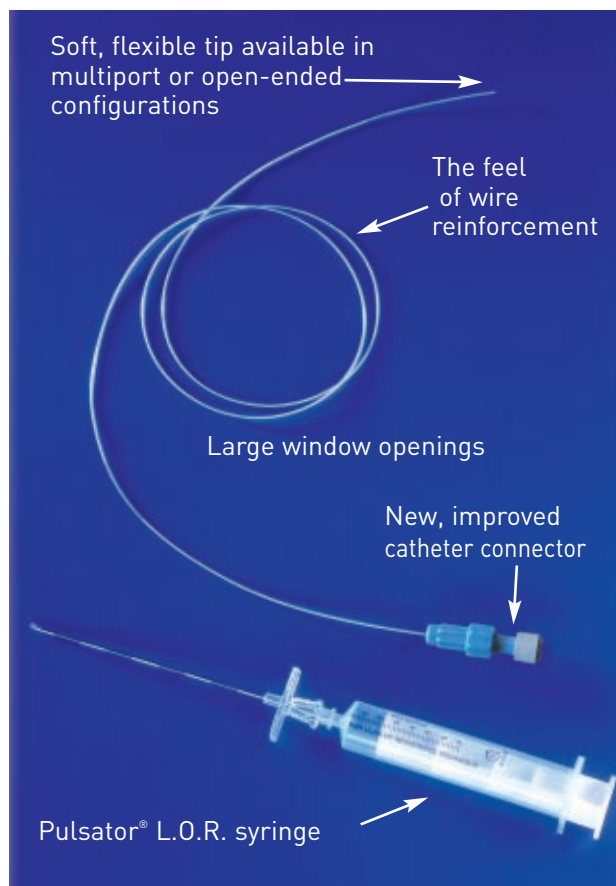

Ordering Information

DuraFlex™ and DuraFlex™ Plus Epidural Catheters

DuraFlex™ Epidural Catheters

REF. CODE	DESCRIPTION	UNITS /CASE
4917-19	DuraFlex™ Open End Epidural Catheter	10
4917M-17	DuraFlex™ Closed End Multiport Epidural Catheter	10

DuraFlex™ Plus Epidural Catheters

REF. CODE	DESCRIPTION	UNITS /CASE
4917S-19	DuraFlex™ Plus Open End Epidural Catheter	10
4917SM-19	DuraFlex™ Plus Closed End Multiport Epidural Catheter	10

All individually packaged, sterile DuraFlex™ catheters are spiral-wound wire reinforced, nylon-coated for use with 16g or 17g epidural needles. Additionally, they include a catheter connector, catheter insertion device, and an epidural catheter label.

Contact your Smiths Medical sales representative today to have a new tray with the DuraFlex™ or DuraFlex™ Plus catheter you choose configured to meet your specific needs.

Customer Service: 1-800-258-5361

www.smiths-medical.com

smiths

Smiths Medical MD, Inc.
Medication Delivery and
Patient Monitoring Division
St. Paul, MN 55112, USA
www.smiths-medical.com

**Rx
ONLY**

Please refer to the Instructions for Use supplied with the product for detailed instructions, warnings, and cautions. DuraFlex, Pulsator and the Smiths design mark are trademarks of the Smiths Medical family of companies. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. ©2006 Smiths Medical family of companies. All rights reserved. 3/06 19349