

## Easy Glide™

sheath introducers for angiographic and interventional procedures

### Features and Benefits

#### Smooth and easy insertion - "Easy Glide"

- balanced proportion of stability and flexibility
- easy atraumatic insertion
- extremely well engineered taper of the sheath onto dilator
- long rounded dilator tip reduces any friction during insertion

#### Smooth catheter loading and control

- smooth guidewire and catheter loading
- easy catheter torquing and manipulation

#### Effective haemostasis


- improved and advanced haemostasis valve design
- prevents any bleed back with the use of catheters

#### Kink protection and resistance


- special flexible strain relief hub sleeve
- sheath shaft material resists permanent kinking

#### "Click" connection

- engaging dilator and sheath
- offers maximum safety

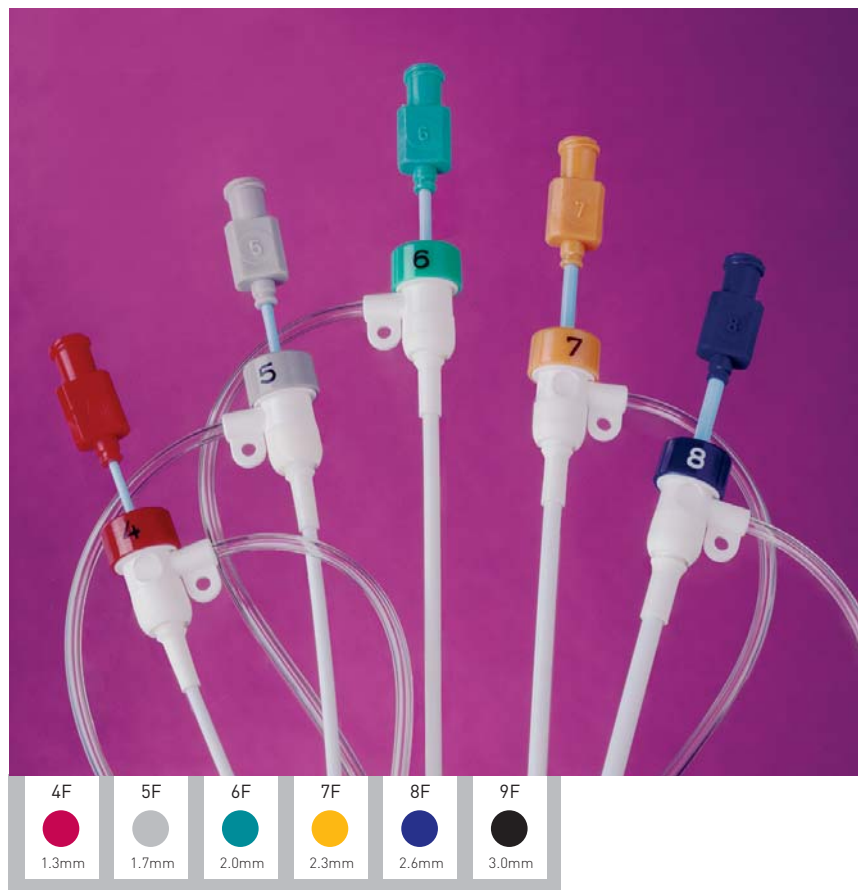

INTERVENTIONALIMAGING

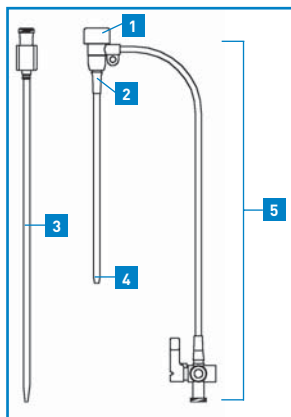

- 1 Sheath**
  - haemostasis valve prevents any bleed back and protects against air embolus
  - catheter manipulation through the valve is simplified
  - suture wing allows secure anchoring
  - colour-coded size indication
- 2 Kink Protection**
  - flexible link prevents the shaft from kinking
- 3 Dilator**
  - facilitates atraumatic catheter insertion
  - aids enlargement of the puncture site without discomfort
  - colour-coded size indication
- 4 Sheath Shaft**
  - removes any ridges at the interface
  - ensures trouble free insertion of the sheath
  - minimises the risk of thrombosis
- 5 Sideport**
  - suitable for both invasive pressure monitoring and for the injection of fluids and taking of fluid samples

## General Specification

	Standard	Long	Paediatric
Sheaths	110mm	230mm	90mm
Dilators	180mm	310mm	140mm
Introducer Needle	18G, 7cm length	18G, 7cm length	20G, 7cm length
Guidewire (4F)	0.89mm [0.035"], 50cm length		
Guidewire (5F - 9F)	0.97mm [0.038"], 50cm length	0.97mm [0.038"], 70cm length	0.64mm [0.025"], 50cm length
Obturator		275mm (120 Series only)	

## Easy Glide™ Kits

Product Code	Description	Qty
<b>100 Series</b> (sheath with haemostasis valve, sideport, dilator, guidewire and insertion needle)		
MXI 100-10-40-35	Introducer Sheath 4F	20
MXI 100-10-50	Introducer Sheath 5F	20
MXI 100-10-60	Introducer Sheath 6F	20
MXI 100-10-70	Introducer Sheath 7F	20
MXI 100-10-80	Introducer Sheath 8F	20
MXI 100-10-90	Introducer Sheath 9F	20
<b>102 Series</b> (sheath with haemostasis valve, sideport, dilator and guidewire)		
MXI102-10-40-35	Introducer Sheath 4F	20
MXI102-10-50	Introducer Sheath 5F	20
MXI102-10-60	Introducer Sheath 6F	20
MXI102-10-70	Introducer Sheath 7F	20
MXI102-10-80	Introducer Sheath 8F	20
MXI102-10-90	Introducer Sheath 9F	20
<b>108 Series</b> (sheath with haemostasis valve, sideport and dilator)		
MXI108-10-40-35	Introducer Sheath 4F	20
MXI108-10-50	Introducer Sheath 5F	20
MXI108-10-60	Introducer Sheath 6F	20
MXI108-10-70	Introducer Sheath 7F	20
MXI108-10-80	Introducer Sheath 8F	20
MXI108-10-90	Introducer Sheath 9F	20
<b>102 Series Long</b> (long sheath with haemostasis valve, sideport, dilator, and guidewire)		
MXI102-24-50	Introducer Sheath 5F	5
MXI102-24-60	Introducer Sheath 6F	5
MXI102-24-70	Introducer Sheath 7F	5
MXI102-24-80	Introducer Sheath 8F	5
MXI102-24-90	Introducer Sheath 9F	5
<b>108 Series Long</b> (long sheath with haemostasis valve, sideport, and dilator)		
MXI108-24-50	Introducer Sheath 5F	5
MXI108-24-60	Introducer Sheath 6F	5
MXI108-24-70	Introducer Sheath 7F	5
MXI108-24-80	Introducer Sheath 8F	5
MXI108-24-90	Introducer Sheath 9F	5
<b>120 Series Long</b> (long sheath with haemostasis valve, sideport, dilator, guidewire, insertion needle, and obturator)		
MXI120-24-50	Introducer Sheath 5F	5
MXI120-24-60	Introducer Sheath 6F	5
MXI120-24-70	Introducer Sheath 7F	5
MXI120-24-80	Introducer Sheath 8F	5
MXI120-24-90	Introducer Sheath 9F	5
<b>Paediatrics</b> (sheath with haemostasis valve, sideport, dilator, guidewire and insertion needle)		
MXI100-09-40	Introducer Sheath 4F	20
MXI102-09-40	Introducer Sheath 4F	20
MXI108-09-40	Introducer Sheath 4F	20

THE DETAILS GIVEN IN THIS LEAFLET ARE CORRECT AT THE TIME OF GOING TO PRESS. THE COMPANY RESERVES THE RIGHT TO IMPROVE THE EQUIPMENT SHOWN.

For further information please call Smiths Medical Customer Services Direct on +44 (0)1923 241411


# smiths

## Smiths Medical International Ltd

Colonial Way, Watford, Herts WD24 4LG, UK  
 Tel: +44 (0)1923 246434 Fax: +44 (0)1923 231595  
[www.smiths-medical.com](http://www.smiths-medical.com)

Easy Glide and the Smiths design mark are trademarks of the Smiths Medical family of companies. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries.

Literature No. TPF 00754-1  
 © 2006 Smiths Medical family of companies. All rights reserved.



# smiths

Smiths Medical - a part of Smiths Group plc