


Portex® Choice Tracheal Tubes for the Operating Room
Affordable quality - SMILE!



AIRWAY MANAGEMENT

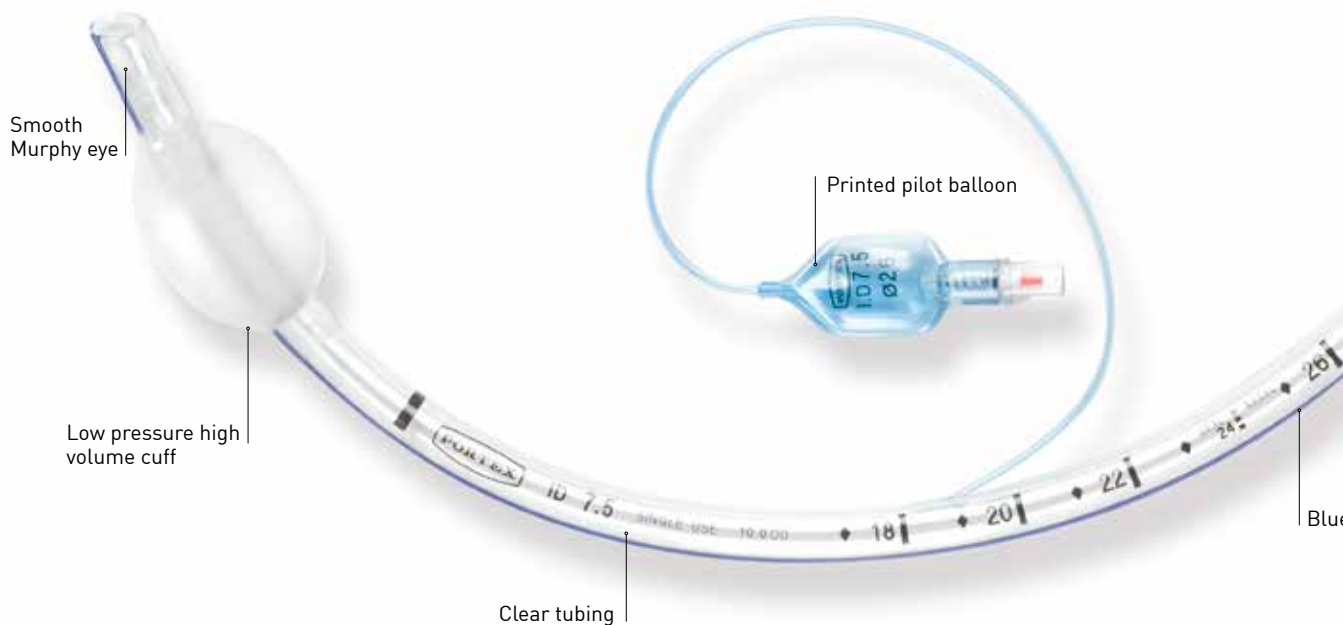

The New Portex® Tracheal Tube

Smiths Medical, a leading global provider of medical devices, is pleased to announce the launch of the new Portex® Tracheal Tube.

The new Portex® Tracheal Tube is the first choice tube for everyday intubation, meeting all the clinical requirements expected of a short-term tube for the operating room.

Affordable quality

Portex®, a renowned brand in airway products, has created a high quality tracheal tube that incorporates all the superior features expected from a Smiths Medical product, yet at an affordable price.

The new Portex® Tracheal Tube is perfect for the majority of all operating room procedures of a short duration.

Low pressure seal

The new Portex® Tracheal Tube has a low pressure, high volume barrel cuff, reducing the risk of creasing on inflation and in turn minimising the pressure on the tracheal wall.

The large cuff resting diameter relieves the pressure on the tracheal wall ensuring patient safety and comfort.



Low pressure high volume cuff



Intubation depth marker



Smooth Murphy eye



Printed pilot balloon

Superior smooth finish

Unlike some competitor tubes that contain a sharp, cold cut edge to the Murphy eye, the new Portex® Tracheal tube contains a smooth and polished Murphy eye.

It also has been smoothed and fully blended to the cuff edge to minimise the risk of trauma to the patient during intubation.

Easily identifiable

A blue radio-opaque line can be seen on the full length of the new Portex® Tracheal Tube to allow the verification of tube positioning within the patient via X-ray.

The tube is made out of clear, implant-tested PVC, therefore correct intubation is confirmed when the tube 'mists'.

All markings on the tube are clearly printed and include an intubation depth marker, vocal cord distal markings and an oral/nasal cut line.

Fully informed

The new Portex® Tracheal tube incorporates a smooth edged pilot balloon that displays the size and maximum inflation pressure allowing for the relevant size to be identified, even when in situ. The balloon changes shape when inflated and contains a central channel so it will not collapse if the cuff needs to be deflated quickly in an emergency situation.

The inflation line is blue to reduce the risk of being accidentally cut.

Portex® Choice Tracheal Tube

For short-term intubation within the operating room

ORDERING INFORMATION				
Product Code	I.D.(mm)	O.D.(mm)	Cuff resting diam.(mm)	Pack Qty
100/150/050	5.0	6.7	22	10
100/150/055	5.5	7.3	22	10
100/150/060	6.0	8.0	24	10
100/150/065	6.5	8.5	24	10
100/150/070	7.0	9.2	26	10
100/150/075	7.5	10.0	26	10
100/150/080	8.0	10.7	28	10
100/150/085	8.5	11.3	28	10
100/150/090	9.0	12.0	30	10
100/150/095	9.5	12.7	30	10



Portex® Soft Seal® Tracheal Tube

For long-term intubation and short-term ICU stays (up to 2 days)

ORDERING INFORMATION				
Product Code	I.D.(mm)	O.D.(mm)	Cuff resting diam.(mm)	Pack Qty
100/199/050	5.0	6.9 mm	25	20
100/199/055	5.5	7.6 mm	28	20
100/199/060	6.0	8.2 mm	29	20
100/199/065	6.5	8.9 mm	30	20
100/199/070	7.0	9.6 mm	31	20
100/199/075	7.5	10.3 mm	32	20
100/199/080	8.0	10.9 mm	33	20
100/199/085	8.5	11.6 mm	33	20
100/199/090	9.0	12.3 mm	33	20
100/199/095	9.5	13.0 mm	33	20



Portex® Suction Above The Cuff Tracheal Tube (SACETT™)

For all ICU patients over 2 days to help reduce the risk of V.A.P

ORDERING INFORMATION				
Product Code	I.D.(mm)	O.D.(mm)	Cuff resting diam.(mm)	Pack Qty
100/189/060	6.0	9.0	23	10
100/189/065	6.5	9.7	23	10
100/189/070	7.0	10.4	30	10
100/189/075	7.5	11.1	30	10
100/189/080	8.0	11.9	30	10
100/189/085	8.5	12.4	30	10
100/189/090	9.0	12.8	30	10



THE DETAILS GIVEN IN THIS LEAFLET ARE CORRECT AT THE TIME OF GOING TO PRESS. THE COMPANY RESERVES THE RIGHT TO IMPROVE THE EQUIPMENT SHOWN

Smiths Medical International Ltd
 1500, Eureka Park, Lower Pemberton, Ashford, Kent TN25 4BF
 Tel: +44 (0)1303 260551 Fax: +44 (0)1303 266761

Smiths Medical ASD, Inc.
 6000 Nathan Lane North, Minneapolis, MN 55442, USA
 Tel: 1-614-210-7300 Toll-Free USA: 1-800-258-5361

www.smiths-medical.com

Soft Seal, SACETT, Portex and the Smiths Medical Design Mark are trademarks of the Smiths Medical family of companies.
 © 2019 Smiths Medical family of companies. All rights reserved.
 Literature No. GA194010GL_122019



smiths medical