

Portex® Blue Line Ultra® Suctionaid®

For easier tracheostomy maintenance



AIRWAY MANAGEMENT

PORTEX®

Blue Line Ultra® Suctionaid®

A new addition to the Blue Line Ultra® range, the Blue Line Ultra Suctionaid® is ideal for improving patient well-being by maintaining a clean, hygienic and unobstructed airway and having the ability to remove secretions from above the cuff.

- Reducing the risk of infection
- Reducing the risk of aspiration
- Maintaining patient comfort
- All the benefits of Blue Line Ultra® including Soft Seal cuff - low pressure, high volume to minimise risk of tracheal trauma while maintaining an effective seal

Reducing the potential risk of infection

Accumulation and stasis of contaminated mucus and subglottic secretions above the cuff of tracheostomy tubes can be uncomfortable for patients and provides an ideal growth medium for pathogens. The integral suction lumen of Blue Line Ultra Suctionaid allows removal of these pooled secretions.

Reducing the risk of aspiration

However good the cuff, microaspiration of contaminated material can potentially lead to pulmonary infection. Maintaining tracheostomy hygiene and regular removal of secretions with Blue Line Ultra Suctionaid can help reduce this aspiration.


Maintaining patient comfort

Blue Line Ultra® Suctionaid® makes it easier to maintain patient comfort.


- Inner Cannula - can be removed, cleaned and re-inserted without disturbing the tracheostomy tube itself
- Tracheostomy tube may be left in situ for up to 30 days reducing the frequency of tube changes
- Suction lumen integrated into tube wall for smooth, easy insertion to minimise trauma
- Tube material softens at body temperature to conform to patient anatomy
- Clear, flexible flange for comfort and aesthetic acceptability




Suction Above the Cuff - Suctionaid[®]

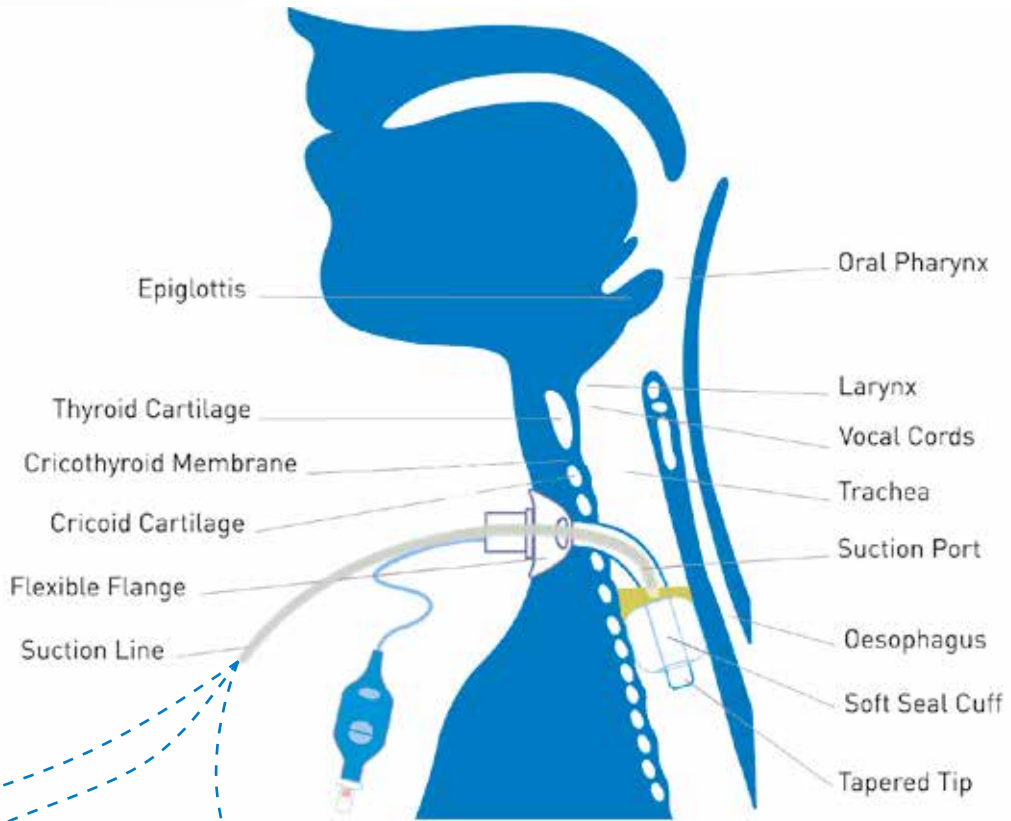


Syringe Aspiration

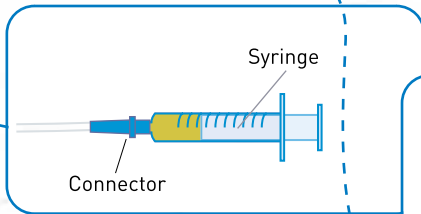

Intermittent Suction

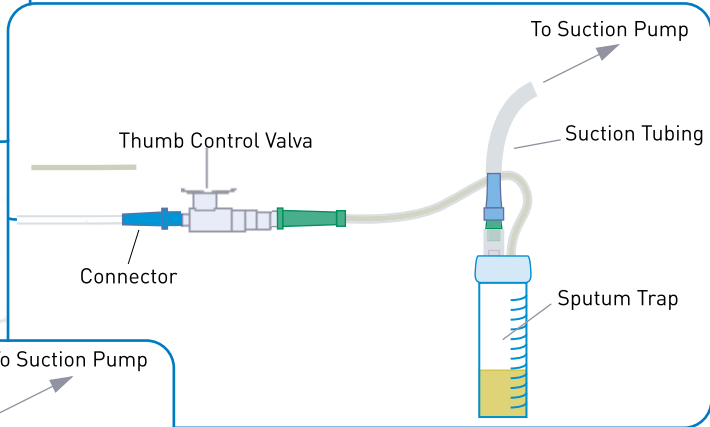

Continuous Suction

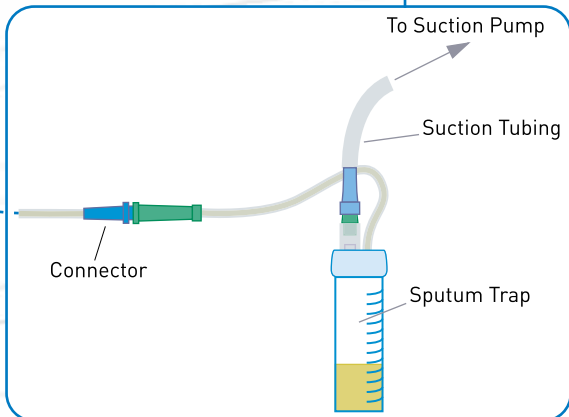

When using continuous or intermittent suctioning, use low level suction up to a maximum of 300mmHg

Ordering Information



100/860

Blue Line Ultra® Suctionaid® Tracheostomy Tube

PRODUCT CODE	TUBE SIZE	TUBE I.D.	TUBE O.D.	TUBE LENGTH	UNITS PER CASE
100/860/060	6.0 mm	6.0 mm	9.2 mm	64.5 mm	10
100/860/070	7.0 mm	7.0 mm	10.5 mm	70.0 mm	10
100/860/075	7.5 mm	7.5 mm	11.3 mm	73.0 mm	10
100/860/080	8.0 mm	8.0 mm	11.9 mm	75.5 mm	10
100/860/085	8.5 mm	8.5 mm	12.6 mm	78.0 mm	10
100/860/090	9.0 mm	9.0 mm	13.3 mm	81.0 mm	10
100/860/100	10.0 mm	10.0 mm	14.0 mm	87.5 mm	10



100/870

Blue Line Ultra® Suctionaid® Tracheostomy Tube Kit

PRODUCT CODE	TUBE SIZE	TUBE I.D.	TUBE O.D.	TUBE LENGTH	UNITS PER CASE
100/870/060	6.0 mm	6.0 mm	9.2 mm	64.5 mm	1
100/870/070	7.0 mm	7.0 mm	10.5 mm	70.0 mm	1
100/870/075	7.5 mm	7.5 mm	11.3 mm	73.0 mm	1
100/870/080	8.0 mm	8.0 mm	11.9 mm	75.5 mm	1
100/870/085	8.5 mm	8.5 mm	12.6 mm	78.0 mm	1
100/870/090	9.0 mm	9.0 mm	13.3 mm	81.0 mm	1
100/870/100	10.0 mm	10.0 mm	14.0 mm	87.5 mm	1

Includes inner cannula with disconnection wedge, tracheostomy tube holder, cleaning brush, and patient label.



100/850

100/851

Inner Cannula for Blue Line Ultra® Tracheostomy Tubes

PRODUCT CODE UNFENESTRATED	PRODUCT CODE FENESTRATED	DESCRIPTION	INNER CANNULA I.D.	UNITS PER CASE
100/850/060	100/851/060	Inner Cannula for 6.0 mm Tube	5.0 mm	2
100/850/070	100/851/070	Inner Cannula for 7.0 mm Tube	6.0 mm	2
100/850/075	100/851/075	Inner Cannula for 7.5 mm Tube	6.5 mm	2
100/850/080	100/851/080	Inner Cannula for 8.0 mm Tube	7.0 mm	2
100/850/085	100/851/085	Inner Cannula for 8.5 mm Tube	7.5 mm	2
100/850/090	100/851/090	Inner Cannula for 9.0 mm Tube	8.0 mm	2
100/850/100	100/851/100	Inner Cannula for 10.0 mm tube	8.5 mm	2



100/891

Blue Line Ultra® Suctionaid® Percutaneous Tracheostomy Kit with Forceps

PRODUCT CODE	TUBE SIZE	TUBE I.D.	TUBE O.D.	TUBE LENGTH	UNITS PER CASE
100/891/070	7.0 mm	7.0 mm	10.5 mm	70.0 mm	1
100/891/080	8.0 mm	8.0 mm	11.9 mm	75.5 mm	1
100/891/090	9.0 mm	9.0 mm	13.3 mm	81.0 mm	1



100/893

Blue Line Ultra® Suctionaid® Percutaneous Tracheostomy Kit without Forceps

PRODUCT CODE	TUBE SIZE	TUBE I.D.	TUBE O.D.	TUBE LENGTH	UNITS PER CASE
100/893/070	7.0 mm	7.0 mm	10.5 mm	70.0 mm	1
100/893/080	8.0 mm	8.0 mm	11.9 mm	75.5 mm	1
100/893/090	9.0 mm	9.0 mm	13.3 mm	81.0 mm	1

Smiths Medical Danmark ApS
 Slotsmarken 18,2970 Hørsholm
 Tel. 7027 2095
 Fax.7027 2098
www.smiths-medical.com

Smiths Medical Sverige AB
 Box 521,192 05 Sollentuna
 Tel. 08- 594 772 50
 Fax.08-594 772 69
www.smiths-medical.com

Smiths Medical ASD, Inc.
 6000 Nathan Lane North
 Minneapolis, MN 55442, USA
 Tel: 1-614-210-7300 Toll-Free USA: 1-800-258-5361
www.smiths-medical.com

Smiths Medical is part of the global technology business Smiths Group plc. Product(s) described may not be licensed or available for sale in all countries. Portex, Blue Line Ultra, Suctionaid and the Smiths Medical and Portex design marks are trademarks of the Smiths Medical family of companies. The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. ©2019 Smiths Medical. All rights reserved. MX194461SEDK-122019



smiths medical