

# smiths medical

## level 1®



### **EQUATOR® Convective Warming System**


Surgery has its complications.  
Take hypothermia off the list.

The EQUATOR® convective warming system is designed to help maintain a patient's normothermic temperature.



# WARMER PATIENTS. BETTER OUTCOMES.

The EQUATOR® convective warming system is designed to maintain patient normothermia and provide thermal comfort before, during and after surgery. Patented hose-end temperature control facilitates consistent temperature of delivered air regardless of ambient temperature. An over-temperature alarm at each set-point (36° C, 40° C and 44° C) promotes patient safety. Four hose-end openings maintain air flow if the hose end is ever occluded by blanket material, helping prevent heat build-up in a single area. A right-angle hose connector also provides gentle entry into the blanket.

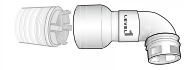

## EQUATOR® Convective Warmer Specifications

Product Number	EQ-5000	SWU-9001	SWU-9002	SWU-9003	SWU-9004
Size	30 cm H x 24 cm W x 19 cm D (11.75" H x 9.5" W x 7.5" D)	For hospitals that have already purchased other manufacturer's convective warmers and still want to enjoy the benefits of Snuggle Warm® convective warming blankets, Snuggle Warm® Hose Adapters are also available.			
Weight	6.8 kg (15 lbs)				
Electrical Configurations	100 VAC, 50/60 Hz 115 VAC, 50/60 Hz 230 VAC, 50 Hz				
Accessories	Air Filter (F3-5000) Sheet Clip (SC-5000) Cord Wrap (CW-5000) Hose (SW5-HOSE-7)				

### References

- Horn EP, et al., Active warming during cesarean delivery. *Anesth Analg.* 2002; 94:409-414.

PRODUCT(S) DESCRIBED MAY NOT BE LICENSED OR AVAILABLE FOR SALE IN CANADA AND OTHER COUNTRIES

**Smiths Medical ASD, Inc.**  
6000 Nathan Lane North  
Minneapolis, MN 55442, USA  
Tel: 1-614-210-7300  
Toll-Free USA: 1-800-258-5361  
[www.smiths-medical.com](http://www.smiths-medical.com)

Find your local contact information at: [www.smiths-medical.com/customer-support](http://www.smiths-medical.com/customer-support)

Smiths Medical is part of the global technology business Smiths Group plc. Please see the Instructions for Use/Operator's Manual for a complete listing of the indications, contraindications, warnings and precautions. Equator, Snuggle Warm, Level 1 and the Smiths Medical design mark are trademarks of Smiths Medical. The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries ©2019 Smiths Medical. All rights reserved. TM194869GB-072019

**Rx ONLY** **CE 2797**

MMSPCA-0230

**smiths medical**