

smiths medical
jelco®



EDGE™ Safety Hypodermic Devices

Safety Systems



LEADING EDGE IN NEEDLE SAFETY

With a comprehensive range of sizes available as stand-alone safety needles and as combinations, EDGE™ needles support standardization across clinical areas and procedures. EDGE™ needle technology also combines the benefits of an advanced safety hypodermic needle with the properties of a conventional needle.

As with conventional hypodermic needles, the stated length of an EDGE™ needle is “true” which translates into no added dead space. Efficient use of EDGE™ needles yields medication savings when compared to other safety hypodermic needles.

EDGE™ COMBINATIONS – SAFETY DEVICE AND SYRINGE ALL IN ONE

1 mL, 3 mL, 5 mL and 10 mL combinations



SINGLE-HAND ACTIVATION KEEPS YOU SAFE

With EDGE™ needles, fingers stay well behind the bevel while activating the safety sheath, helping reduce the potential for needlestick injuries by minimizing exposure. With other devices, the greatest exposure to possible needle-stick injuries occurs when activating the safety mechanism.




Rest the tip of the safety sheath on any flat surface with the needle positioned above.

INTEGRAL SAFETY WITH THE LOOK OF A CONVENTIONAL NEEDLE

Molded arrow points in
direction of bevel

ISO standard
colour-coded hub
and safety sheath
assembly enables
rapid identification of
needle gauge size

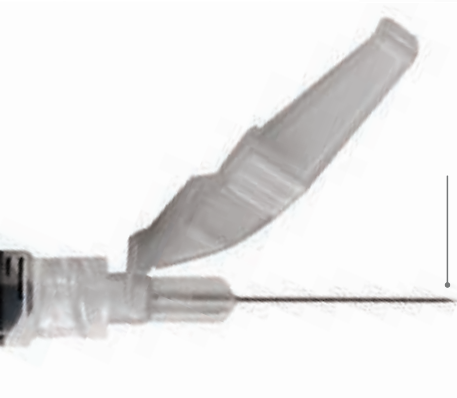

Bevel-oriented needle
– bevel is “up” when
the safety sheath is
positioned at 3-o’clock
so the injection site is
not obscured



Apply downward pressure on
the syringe until the needle
locks into the safety sheath
with an audible click.



Discard the syringe and needle
in accordance with sharp safety
procedures.

EDGE™ Safety Device

Code	Description	Units Per Case
401810	18g x 1"	1000
401815	18g x 1 ½"	1000
401910	19g x 1"	1000
401915	19g x 1 ½"	1000
402010	20g x 1"	1000
402015	20g x 1 ½"	1000
402110	21g x 1"	1000
402115	21g x 1 ½"	1000
402210	22g x 1"	1000
402215	22g x 1 ½"	1000
402310	23g x 1"	1000
402315	23g x 1 ½"	1000
402558	25g x 5/8"	1000
402510	25g x 1"	1000
402515	25g x 1 ½"	1000
402605	26g x ½"	1000
402705	27g x ½"	1000
4027125	27g x 1 ¼"	1000
403005	30g x ½"	1000

EDGE™ Safety Combinations*

431810	18g x 1", 3 mL syringe	400
431815	18g x 1 ½", 3 mL syringe	400
431910	19g x 1", 3 mL syringe	400
431915	19g x 1 ½", 3 mL syringe	400
432010	20g x 1", 3 mL syringe	400
432015	20g x 1 ½", 3 mL syringe	400
432110	21g x 1", 3 mL syringe	400
432115	21g x 1 ½", 3 mL syringe	400
432210	22g x 1", 3 mL syringe	400
432215	22g x 1 ½", 3 mL syringe	400
432310	23g x 1", 3 mL syringe	400
432558	25g x 5/8", 3 mL syringe	400
432510	25g x 1", 3 mL syringe	400
452010	20g x 1", 5 mL syringe	400
452015	20g x 1 ½", 5 mL syringe	400
452110	21g x 1", 5 mL syringe	400
452115	21g x 1 ½", 5 mL syringe	400
452215	22g x 1 ½", 5 mL syringe	400
4102010	20g x 1", 10 mL syringe	400
4102015	20g x 1 ½", 10 mL syringe	400
4102115	21g x 1 ½", 10 mL syringe	400

*Luer Lock Syringe

EDGE™ INSULIN AND TB INJECTIONS

EDGE™ Permanently Attached Insulin Syringes*

Code	Description	Units Per Case
4428-5	0.5 mL, 28g x ½"	600
4428-1	1 mL, 28g x ½"	600
4429-3	0.3 mL, 29g x ½"	600
4429-5	0.5 mL, 29g x ½"	600
4429-1	1 mL, 29g x ½"	600
4430-3	0.3 mL, 30g x ½"	600
4430-5	0.5 mL, 30g x ½"	600

* For U-100 Insulin

EDGE™ Detachable Insulin Syringes*

4317	1 mL, 26g x ½"	400
4319	1 mL, 27g x ½"	400

* Luer Slip Needle For U-100 Insulin

EDGE™ Permanently Attached TB Syringes

4425-TB	1 mL, 25g x 5/8"	600
4426-TB	1 mL, 26g x 3/8" Intradermal Bevel	600
4427-TB	1 mL, 27g x ½"	600

EDGE™ Detachable TB Syringes†

412558	25g x 5/8"	400
412705	27g x ½"	400

†Luer Lock Needle

Blunt Fill Needle

BN1815	18g x 1½"	1000
--------	-----------	------



Customer Support (0)845 850 0445
www.smiths-medical.com

PRODUCT(S) DESCRIBED MAY NOT BE LICENSED OR AVAILABLE FOR SALE IN CANADA AND OTHER COUNTRIES

Smiths Medical International Ltd.

1500 Eureka Park, Lower Pemberton
Ashford, Kent, TN25 4BF
Tel: +44 (0)845 850 0445

MKEECA-0167

smiths medical

www.smiths-medical.com

Find your local contact information at
<http://www.smiths-medical.com/customer-support/contact-us>

Smiths Medical is part of the global technology business Smiths Group plc. Please see the Instructions for Use/Operator's Manual for a complete listing of the indications, contraindications, warnings and precautions. Jelco and the Smiths Medical design mark are trademarks of Smiths Medical. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries.

©2017 Smiths Medical. All rights reserved. SS195202GB-112017