


smiths medical
deltec®

**PORT-A-CATH® and P.A.S. PORT® T2 POWER P.A.C.® Systems
GRIPPER PLUS® POWER P.A.C. Safety Huber Needles**

The proven performance you expect, now in a powerful system

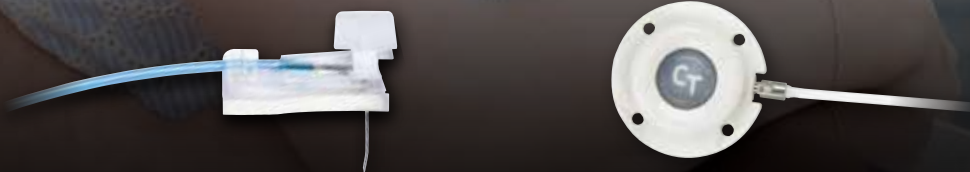

A POWERFUL SYSTEM WITH POWERFUL FEATURES

We have the largest variety and most complete system of ports and needles for power injection for any therapy need allowing you to choose the optimal solution for your patients' needs.

A COMPLETE OFFERING OF POWER-INJECTABLE PORT SYSTEMS

Choice of:

Materials – polysulfone or titanium

Profile – standard or low profile

Lumen configuration – single or dual

Catheter options – multiple French sizes

Introducer – regular peel-away and valved peel-away

PATIENT IDENTIFICATION SYSTEM

Patients are provided with identification and educational materials for their PORT-A-CATH® and P.A.S. PORT® T2 POWER P.A.C.® systems.

- Identification cards
- Key ring identifier
- Patient information guide
- Medical record stickers





29 different portal configurations

PORT-A-CATH® AND P.A.S. PORT® T2 POWER P.A.C.® SYSTEMS

Proven technology for over 25 years:

- Radiopaque CT embedded into port septum and able to withstand thousands of needlesticks while remaining visible under x-ray
- No change in power injection performance after hundreds of needle punctures
- Allows power injections up to 5 mL/sec and 300 psi pressure rating
- MR-conditional up to 3 Tesla
- Compatible with all power needles
- Less thrombosis than other ports¹
- Low complication rates²
- Superior PolyFlow® polyurethane catheter

14 different needle configurations

GRIPPER PLUS® POWER P.A.C. NEEDLES

Designed for power injections by allowing higher pressure and flow rates with power-injectable port systems.

Innovative features have set the standard in reliable portal access:

- Blue tubing and clamp markings for easy power injection identification
- Allows power injections up to 5 mL/sec and 300 psi pressure rating
- Compatible with all power ports
- Unique safety arm provides audible “click” to ensure needle is secured
- Non-coring and non-siliconized needle
- Nonabsorbent, closed foam pad for increased patient comfort
- Fully threaded Luer-lock connector fits standard accessories

GRIPPER PLUS® POWER P.A.C. SAFETY HUBER NEEDLES PRODUCT SPECIFICATIONS

Gauge	Length	Maximum Flow Rate	Reorder # Nac	Reorder # Non Y-Site
19	3/4 inch	5 ml/sec	21-3468-24	21-3368-24
19	1 inch	5 ml/sec	21-3469-24	21-3364-24
19	1 1/4 inch	5 ml/sec	21-3470-24	21-3365-24
20	5/8 inch	5 ml/sec	21-3464-24	21-3369-24
20	3/4 inch	5 ml/sec	21-3465-24	21-3367-24
20	1 inch	5 ml/sec	21-3466-24	21-3362-24
20	1 1/4 inch	5 ml/sec	21-3467-24	21-3363-24

REFERENCES

1. Strategic Healthcare Programs, LLC. Catheter-related complications by catheter type. 2004.
2. Biffi, Robert M.D. Use of totally implantable central venous access ports for high-dose chemotherapy and peripheral blood stem cell transplantation: results of a monocentre series of 376 patients. *Annals of Oncology* 2004; 15:296-300.

WE HAVE THE LARGEST VARIETY AND MOST COMPLETE SYSTEM OF PORTS AND NEEDLES FOR POWER INJECTION FOR ANY THERAPY NEED.

PORT-A-CATH® AND P.A.S. PORT® T2 POWER P.A.C.® SYSTEMS PRODUCT SPECIFICATIONS

PORTAL

CATHETER

Single-Lumen Standard

	Material	Base (mm)	Height (mm)	Weight (g)	Septum Diameter (mm)	Connector	Material	ID (mm)	OD (mm)	French (Fr)	Length (cm)	Introducer Size (Fr)	Configuration	Reorder #
PORT-A-CATH® POWER P.A.C.	Titanium	30.5	14.7	16	11.4	Ultra-Lock®	PolyFlow® Polyurethane	1.6	2.6	7.8	76	8.5	Unassembled	Kit: 21-4422-24 Tray: 21-4423-24
PORT-A-CATH® POWER P.A.C.	Titanium	30.5	14.7	16	11.4	Ultra-Lock®	PolyFlow® Polyurethane	1.0	1.9	5.8	76	6	Unassembled	Kit: 21-4424-24 Tray: 21-4425-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	28	14.6	16	11.4	Ultra-Lock®	PolyFlow® Polyurethane	1.0	1.9	5.8	76	6	Unassembled	Kit: 21-4452-24 Tray: 21-4453-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	28	14.6	16	11.4	Ultra-Lock®	PolyFlow® Polyurethane	1.6	2.6	7.8	76	8.5	Unassembled	Kit: 21-4454-24 Tray: 21-4455-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	28	14.6	16	11.4	Ultra-Lock®	PolyFlow® Polyurethane	1.0	1.9	5.8	76	6	Assembled	Tray: 21-4463-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	28	14.6	16	11.4	Ultra-Lock®	PolyFlow® Polyurethane	1.6	2.6	7.8	76	8.5	Assembled	Tray: 21-4465-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	30.5	14.7	10	11.4	Ultra-Lock®	PolyFlow® Polyurethane	1.0	1.9	5.8	76	6*	Unassembled	Tray: 21-4476-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	30.5	14.7	10	11.4	Ultra-Lock®	PolyFlow® Polyurethane	1.6	2.6	7.8	76	8.5*	Unassembled	Tray: 21-4474-24

Single-Lumen Low Profile™

PORT-A-CATH® POWER P.A.C.	Titanium	25	11.5	9.5	9.5	Ultra-Lock®	PolyFlow® Polyurethane	1.0	1.9	5.8	76	6	Unassembled	Kit: 21-4436-24 Tray: 21-4437-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	25	11.5	4.8	9.5	Ultra-Lock®	PolyFlow® Polyurethane	1.0	1.9	5.8	76	6	Unassembled	Kit: 21-4482-24 Tray: 21-4483-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	25	11.5	4.8	9.5	Ultra-Lock®	PolyFlow® Polyurethane	1.0	1.9	5.8	76	6	Assembled	Tray: 21-4485-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	25	11.5	4.8	9.5	Ultra-Lock®	PolyFlow® Polyurethane	1.6	2.6	7.8	76	8.5	Unassembled	Kit: 21-4470-24 Tray: 21-4471-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	25	11.5	4.8	9.5	Ultra-Lock®	PolyFlow® Polyurethane	1.6	2.6	7.8	76	8.5	Assembled	Tray: 21-4473-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	25	11.4	5	9.5	Ultra-Lock®	PolyFlow® Polyurethane	1.0	1.9	5.8	76	6*	Unassembled	Tray: 21-4477-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	25	11.4	5	9.5	Ultra-Lock®	PolyFlow® Polyurethane	1.6	2.6	7.8	76	8.5*	Unassembled	Tray: 21-4475-24

Dual-Lumen Low Profile™

PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	38.7 x23.5	11	11	9.5	Wing-Lock®	PolyFlow® Polyurethane	1.4	3.2	9.6	76	10	Unassembled	Tray: 21-8468-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	38.7 x23.5	11	11	9.5	Wing-Lock®	PolyFlow® Polyurethane	1.0	2.2	6.7	76	7	Unassembled	Tray: 21-8466-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	38.7 x23.5	11	12	9.5	Wing-Lock®	PolyFlow® Polyurethane	1.4/1.4	3.2	9.6	76	10*	Unassembled	Tray: 21-8470-24
PORT-A-CATH® II POWER P.A.C.	Polysulfone/ Titanium	38.7 x23.5	11	12	9.5	Wing-Lock®	PolyFlow® Polyurethane	1.0/1.0	2.2	6.7	76	7*	Unassembled	Tray: 21-8469-24

Peripheral Single-Lumen

P.A.S. PORT® T2 POWER P.A.C.®	Titanium	24.5 x18.2	11.5	8.3	9.5	Ultra-Lock®	PolyFlow® Polyurethane	1.0	1.9	5.8	76	6	Unassembled	Kit: 21-4872-24 Tray: 21-4873-24
-------------------------------	----------	---------------	------	-----	-----	-------------	------------------------	-----	-----	-----	----	---	-------------	-------------------------------------

1. Trays include: portal, catheter, Port-A-Cath® straight needle, blunt needle (except preassembled systems), vein pick, 18ga extra thin wall introducer needle, "J" guidewire(s), vessel dilator/sheath introducer assembly, 90° Port-A-Cath® needle, Gripper Plus® Power P.A.C. safety needle(s), 2 syringes, tunneling tool, and Point-Lok® device.

2. Kits include: portal, catheter, Port-A-Cath® straight needle, blunt needle, vein pick, and Point-Lok® device

* Valved Peel-Away Introducer

For more information visit our website at www.smiths-medical.com

PRODUCT(S) DESCRIBED MAY NOT BE LICENSED OR AVAILABLE FOR SALE IN CANADA AND OTHER COUNTRIES

Smiths Medical ASD, Inc.
6000 Nathan Lane North
Minneapolis, MN 55442, USA
Tel: 1-614-210-7300
Toll-Free USA: 1-800-258-5361
www.smiths-medical.com

MPAUC-1498

smiths medical

Find your local contact information at: www.smiths-medical.com/customer-support

Smiths Medical is part of the global technology business Smiths Group plc. Please see the Instructions for Use/Operator's Manual for a complete listing of the indications, contraindications, warnings and precautions. Port-A-Cath, P.A.S. Port, Low Profile, Ultra-Lock, PolyFlow, Wing-Lock, Point-Lok, Gripper Plus and the Deltec and Smiths Medical design mark are trademarks of Smiths Medical. The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. ©2017 Smiths Medical. All rights reserved. VA192000EN-102017