

SurgiVet® HOTLINE® Blood and Fluid Warmer



The SurgiVet® HOTLINE®

Blood and Fluid Warmer is one of the leading choices in veterinary temperature management. A recirculating solution keeps blood or fluid warm all the way to the patient.

HOTLINE® is fast, requires limited maintenance, and gives you the confidence to place your focus where it's needed—your patient. By keeping patients warm, they are less likely to suffer post-operative shivering, delayed wound healing, delayed recovery from anesthesia and cardiac instability.

For more information, visit www.surgivet.com.

VETERINARY


SurgiVet™



V707000

SurgiVet® HOTLINE® Blood and Fluid Warmer

- The patented patient line with triple lumen tubing maintains a layer of 107.6° F circulating solution around IV line
- Keeps blood and fluids warm between 95-105° F
- Reduces patient line cool down
- Dedicated temperature check device to calibrate the temperature set point
- Minimal set up time and is easy to disconnect from a patient
- Accommodates a slow drip, up to a bolus
- Warm and ready to use in four minutes
- Compact size and lightweight
- Can be mounted to an IV pole
- Remanufactured unit with 1 year warranty
- Minimal maintenance and low operating costs
- Each warmer includes five warming sets
- A leading manufacturer in hypothermia prevention

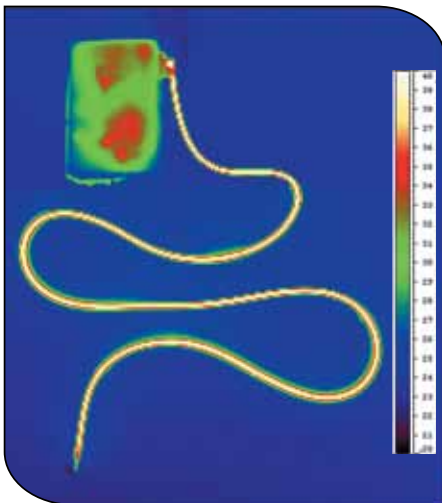

By using thermal photography, you can see a consistent temperature throughout the sample line. No retouching or digital manipulation was used.

Technical Specifications

OPERATING

- **Source Temperature**
106.7° F ± 32° F [41.5° C ± .05° C]
- **Maximum Mounting Height**
42 inches (on a pole)

DIMENSIONS

- 9.5" H x 8.3" W x 7" D
(24 cm H x 21 cm W x 17.7 cm)

WEIGHT

- **Dry**
7.6 lbs. (3.5 kg)
- **Wet**
11 lbs. (5 kg)

ENVIRONMENTAL

- **Operating Temperature**
50° F to 113° F [10° C to 45° C]
- **Storage Temperature**
-4° F to 158° F [-20° C to 70° C]
- **Transport Temperature**
-4° F to 158° F [-20° C to 70° C]

CATALOG

- V707000

Catalog #	Description	Units
V707000	SurgiVet® HOTLINE® Blood & Fluid Warmer	1
32491B9	Blood & Fluid Warming Set with Injection Port	1

SPECIFICATIONS SUBJECT TO CHANGE WITHOUT NOTICE.

Please read the Instructions for Use supplied with the product for detailed instructions, warnings and cautions.

Smiths Medical PM, Inc.

Veterinary
N7W22025 Johnson Drive
Waukesha, Wisconsin 53186 USA
Phone: 262-513-8500 Fax: 262-542-0718
Toll-Free USA: 888-745-6562
www.surgivet.com

Smiths Medical, part of the global technology business Smiths Group

Product(s) described may not be licensed or available for sale outside of the United States.

Hotline, SurgiVet and the Smiths Medical design mark are trademarks of the Smiths Medical family of companies. The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are the trade names, trademarks or service marks of their respective owners. © 2009 Smiths Medical family of companies. All rights reserved. V707000 Rev.00 08/09

Smiths Medical Canada Ltd.

Markham, Ontario, Canada, L3R 4Y8
Phone: 905-477-2000
Toll-Free: 800-387-4346

smiths medical